


Preparing and sending samples

It's essential to use cultures that are as fresh as possible, and that have not been chilled.


Microorganisms on agar medium

- Seal plates with Parafilm
- Label plates/mark colonies for identification
- Wrap it up and send it to us*

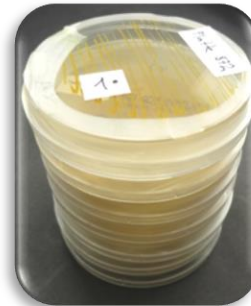

Bacteria and yeasts suspended in 70% ethanol

- Use tubes with screw caps
- Take 1-2 coils from a bacterial colony
- Suspend in 100-200 µl of 70% ethanol (pellet visible)
- Label the tubes
- Wrap it up and send it to us*

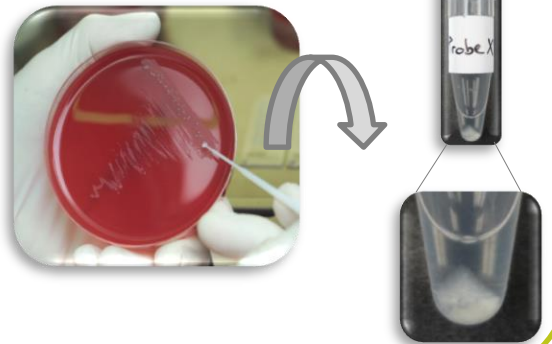

Molds

- Use tubes with screw caps
- Remove 1cm² sample material
- Put it in a tube and label it
- Wrap it up and send it to us*



* Transport regulations: UN3373 logo on the package, as well as the words <<Biological material, CATEGORY B >>